



Xpress™ Minimally Invasive Pedicle Screw System

GENERAL INFORMATION

The Xpress Minimally Invasive Pedicle Screw System consists of rods, pedicle screws, screw caps and hand instruments. Various forms and sizes of these implants are available so that adaptations can be made to take into account the pathology and anatomy of an individual patient. The system implant components are made of Ti6Al4V ELI, a titanium-based alloy which complies with ASTM F136, and Co28Cr6Mo, a cobalt chromium alloy which complies with ASTM F1537.

INDICATIONS FOR USE

The Xpress Minimally Invasive Pedicle Screw System is indicated for the treatment of severe spondylolisthesis (Grade 3 and 4) of the L5-S1 vertebra in skeletally mature patients receiving fusion by autogenous bone graft having implants attached to the lumbar and sacral spine (L3 to sacrum) with removal of the implants after the attainment of a solid fusion.

In addition, the Xpress Minimally Invasive Pedicle Screw System is intended for posterior, non-cervical (T1-S2/ilium) pedicle spinal fixation to provide immobilization and stabilization of the spinal segments in skeletally mature patients as an adjunct to fusion in the treatment of the following acute and chronic instabilities or deformities of the thoracic, lumbar and sacral spine: degenerative spondylolisthesis with objective evidence of neurological impairment, fracture, dislocation, scoliosis, kyphosis, spinal tumor, and failed previous fusion (pseudarthrosis).

CONTRAINDICATIONS

1. Active infectious process in the patient, particularly in or adjacent to the spine or spinal structures.
2. Morbid obesity.
3. Pregnancy.
4. Grossly distorted anatomy due to congenital abnormalities.
5. Any medical or surgical condition which would preclude the potential benefit of spinal implant surgery.
6. Rapid joint disease, bone absorption, osteopenia, osteomalacia, or osteoporosis. Osteopenia or osteoporosis is a relative contraindication since this condition may limit the degree of obtainable correction and/or the amount of mechanical fixation.
7. Suspected or documented metal allergy or intolerance.
8. Any patient in which implant utilization would interfere with anatomical structures or expected physiological performance.
9. Use of these implants is relatively contraindicated in patients whose activity, mental capacity, mental illness, alcohol or drug abuse, occupation or life-style may interfere with their ability to follow post-operative instructions.
10. Any case not needing a bone graft and fusion or where fracture healing is not required.



Xpress™ Minimally Invasive Pedicle Screw System

11. Other medical or surgical conditions which would preclude the potential benefit of spinal implant surgery such as the presence of tumors, congenital abnormalities, elevation of white blood cell count (WBC), or a left shift in the WBC differential count.

⚠ WARNINGS:

1. System implants and instruments are provided non-sterile and must be sterilized before use. Validated sterilization cycle parameter protocols are noted in the STERILIZATION section of this insert.
2. As with all orthopedic and neurosurgical implants, none of the system implants should ever be reused under any circumstances. Risks associated with reuse include infection, non-union (pseudarthrosis), serious patient injury or death.
3. The Xpress Minimally Invasive Pedicle Screw System has not been evaluated for safety and compatibility in the MR environment. It has not been tested for heating, migration, or image artifact in the MR environment. The safety of the Xpress Minimally Invasive Pedicle Screw System in the MR environment is unknown. Scanning a patient who has this device may result in patient injury.
4. The safety and effectiveness of pedicle screw spinal systems have been established only for spinal conditions with significant mechanical instability or deformity requiring fusion with instrumentation. These conditions are significant mechanical instability or deformity of the thoracic, lumbar and sacral spine secondary to severe spondylolisthesis (Grade 3 and 4) of L5-S1 vertebra, degenerative spondylolisthesis with objective evidence of neurological impairment, fracture, dislocation, scoliosis, kyphosis, spinal tumor, and failed previous fusion (pseudoarthrosis). The safety and effectiveness of these devices for any other conditions are unknown.

PRECAUTIONS

The implants must be implanted only by experienced surgeons with specific training in the use of this pedicle screw spinal system because this is a technically demanding procedure presenting a risk of serious injury to the patient. The use of implants must be decided upon with regard to the surgical and medical indications, the potential risks, and limitations related to this type of surgery. The surgeon and patient should demonstrate knowledge of the contraindications, side effects, precautions, metallurgic and biological characteristics of the implants to be used. Xpress implants must not be used together with implants from a different source, a different manufacturer or made from a different material.

Knowledge of surgical techniques, proper reduction, selection and placement of implants, and pre and post-operative patient management are considerations essential to a successful surgical outcome. Appropriate selection, placement and fixation of the spinal system components are critical factors which affect implant service life. As in the case of all prosthetic implants, the durability of these components is affected by numerous biologic, biomechanics and other extrinsic factors, which limit their service life. Accordingly, strict adherence to the indications, contraindications, precautions, and warnings for this



Xpress™ Minimally Invasive Pedicle Screw System

product is essential to potentially maximize service life. (Note: While proper implant selection can minimize risks, the size and shape of human bones present limitations on the size, shape, and strength of the implants).

Patients who smoke have been shown to have an increased incidence of pseudoarthrosis. Such patients should be advised of this fact and warned of the potential consequences.

If the patient is involved in an occupation or activity which applies inordinate stress upon the implant (e.g. substantial walking, running, lifting, or muscle strain) resultant forces can cause failure of the device. In some cases, progression of degenerative disease may be so advanced at the time of implantation that the expected useful life of the appliance may be substantially decreased. In such cases, orthopedic devices may be considered only as a delaying technique or to provide temporary relief. Patients should be instructed in detail about the limitations of the implants, including, but not limited to, the impact of excessive loading through patient weight or activity, and be taught to govern their activities accordingly. The patient should understand that a metallic implant is not as strong as normal, healthy bone and will bend, loosen or fracture if excessive demands are placed on it. An active, debilitated, or demented patient who cannot properly use weight supporting devices may be particularly at risk during postoperative rehabilitation.

Care must be taken to protect the components from being marred, nicked or notched as a result of contact with metal or abrasive objects. Alterations will produce defects in surface finish and internal stresses which may become the focal point for eventual breakage of the implant.

Physician Note: The physician is the learned intermediary between the company and the patient. The indications, contraindications, warnings, and precautions given in this document must be conveyed to the patient.

Benefits of spinal fusion utilizing any pedicle screw fixation system has not been adequately established in patients with stable spine. Potential risks associated with the use of this system, which may require additional surgery, include; device component neurological injury, and vascular or visceral injury. Discard all damaged or mishandled implants. Never reuse an implant, even though it may appear undamaged.

Internal fixation devices cannot withstand activity and load levels equal to those placed on normal healthy bone. Until maturation of the fusion mass is confirmed, do not subject this device to the stress of full weight bearing, or implant failure may result.

Contouring and bending of a system component may reduce its fatigue strength and cause failure under load. If spinal screws are bent or otherwise damaged during insertion or adjustment, they must not be implanted and must be replaced. Rods should only be contoured with the proper contouring instruments.



Xpress™ Minimally Invasive Pedicle Screw System

Incorrectly contoured rods or rods which have been repeatedly or excessively contoured must not be implanted.

Removal of an unloosened spinal screw may require the use of special instruments to disrupt the interface at the implant surface. This technique may require practice in the laboratory before attempted clinically. Any decision by a surgeon to remove the internal fixation device should take into consideration such factors as the risk to the patient of the additional surgical procedure as well as the difficulty of removal. Implant removal should be followed by adequate postoperative management to avoid fracture.

INSTRUCTIONS FOR USE

A successful result is not always achieved in every surgical case, especially in spinal surgery where many extenuating circumstances may compromise results. Preoperative planning and operating procedures, including knowledge of surgical techniques, proper reduction, and proper selection and placement of the implant are critical considerations in achieving a successful result.

For additional instructions regarding the proper use and application of all Xpress implants and instruments, please refer to the Xpress Surgical Technique Manual (available at no charge upon request).

Use of the Xtant Medical Xpress Minimally Invasive Pedicle Screw System should only be considered when the following preoperative, intraoperative and postoperative conditions exist:

Preoperative Management:

1. Only patients that meet the criteria described in the indications should be selected.
2. Patient conditions and/or predispositions such as those mentioned in the contraindications should be avoided.
3. Care should be used in the handling and storage of the implant components. The implants should not be scratched or otherwise damaged. Implants and instruments should be protected during storage, especially from corrosive environments.
4. The type of construct to be assembled for the case should be determined prior to beginning the surgery. An adequate inventory of sizes should be available at the time of surgery, including sizes larger and smaller than those expected to be used.
5. Since mechanical parts are involved, the surgeon should be familiar with the various components before using the equipment and should personally assemble the devices to verify that all parts and necessary instruments are present before the surgery begins. The Xtant Medical Xpress System components are not to be combined with the components from another manufacturer. Different metal types should not be used together.
6. All components and instruments should be cleaned and sterilized prior to use. Additional sterile components should be available in case of unexpected need.



Xpress™ Minimally Invasive Pedicle Screw System

Intraoperative Management:

1. Extreme caution should be used around the spinal cord and nerve roots. Damage to these structures will cause loss of neurological function.
2. Breakage, slippage, or misuse of instruments or implant components may cause injury to the patient or operative personnel.
3. Whenever possible or necessary, an imaging system should be utilized to facilitate surgery.
4. Caution should be taken in handling the implants; Damage to the implants may affect their performance.
5. Implants should not be reused under any circumstances.
6. If the Xtant Medical Xpress system implants are being used in conjunction with Xtant Medical Fortex system components, care must be taken to use the persuader instrument with the gray-colored handle only and not the blue-colored handle.

Postoperative Management:

The physician's postoperative directions and warnings to the patient and the corresponding patient compliance are extremely important.

1. Detailed instructions on the use and limitations of the device should be given to the patient. The patient should be instructed to limit and restrict physical activities, especially lifting and twisting motions and any type of sport participation. Patients should be advised of their inability to bend at the point of spinal fusion and taught to compensate for this permanent restriction in body motion. The patient should be advised not to smoke or consume alcohol during the bone graft healing process.
2. If a non-union develops or the components loosen, bend, and/or break, the device(s) should be revised and/or removed immediately before serious injury occurs. Failure to immobilize a delayed or nonunion of bone will result in excessive and repeated stresses on the implant. By the mechanism of fatigue, these stresses can cause eventual bending, loosening or breakage of the device(s).
3. If required, the device may be disassembled for explantation. Care should be taken to avoid damaging the implant and surrounding tissue as little as possible. The explanted device should be cleaned and disinfected using the instructions provided for cleaning/disinfection of instruments. Information on the procedure and patient should be retained to assist in any investigation.
4. Any retrieved devices should be treated in such a manner that reuse in another surgical procedure is not possible. As with all orthopedic implants, none of the Xpress components should ever be reused under any circumstances.

POTENTIAL COMPLICATIONS AND ADVERSE EFFECTS

1. Early or late loosening of the components.



Xpress™ Minimally Invasive Pedicle Screw System

2. Disassembly, bending or breakage of any or all of the components.
3. Foreign body (allergic) reaction to the implants.
4. Infection.
5. Non-union (pseudarthrosis).
6. Loss of neurological function, including paralysis (complete or incomplete), radiculopathy, dysesthesia, hyperesthesia, anesthesia, paresthesia, development or continuation of pain, numbness, neuroma, tingling sensation, dural tears, neuropathy, neurological deficits (transient, permanent, or delayed), reflex deficits, bilateral paraplegia, and/or arachnoiditis.
7. Hemorrhage, hematoma, seroma, embolism, edema, stroke, excessive bleeding, phlebitis, wound necrosis, or wound dehiscence.
8. Misalignment of anatomical structures or loss of spinal mobility.
9. Bone graft donor complications including pain, fracture or wound healing problems.
10. Atelectasis.
11. Cessation of any potential growth of the operated portion of the spine.
12. Vascular damage resulting in excessive bleeding.
13. Loss or impairment of bowel, sexual, and/or bladder function and other types of urological compromise.
14. Fracture, damage, degenerative changes or instability of any bone above and/or below the level of surgery.
15. Gastrointestinal system compromise.
16. Bone loss due to resorption or stress shielding.
17. Death.

PACKAGING, LABELING AND STORAGE

The implants are supplied clean and NON-STERILE. They must be sterilized (see below). The implants are delivered in packages. These must be intact at the time of receipt. The implants may be delivered as a complete set. Implants and instruments are contained within specially designed trays which can be sterilized directly. Use care in handling and storage of the implant components. Cutting, sharply bending, or scratching the surface can significantly reduce the strength and fatigue resistance of the implant system. This, in turn, could induce cracks and/or non-visible internal stresses that could lead to fracture of the implants. Implants and instruments in storage should be protected from corrosive environments such as salt, air, moisture, etc. Inspection and trial assembly are recommended prior to surgery to determine if instruments or implants have been damaged during the storage processes.

STERILIZATION

All Xpress Minimally Invasive Pedicle Screw System implants and instruments are provided clean and non-sterile and must be sterilized before use. All implants and instruments must be free of packaging material and bio-contaminants prior to sterilization. To achieve a sterility assurance level of not less than 10^{-6} , all non-sterile implants and instruments should be autoclave sterilized using the following validated cycle



Xpress™ Minimally Invasive Pedicle Screw System

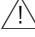
parameters:

Saturated steam method (pressure 41.9 psia, 205.5 kPa), pre-vacuum air removal (3 pulses), 270° F (132° C), 4-minute exposure time, 30-minute drying time, in a double-wrapped case configuration. Note: Most sterilizers automatically control pressure to ensure saturated steam and may not be adjusted by the user.

Use FDA cleared sterilizers and sterilization packaging materials/accessories (wraps, biological indicators, chemical indicators, etc.) for both non-sterile implant components and instruments.

Do Not stack instrument trays during sterilization.

CLEANING OF INSTRUMENTS

 Caution: Use of sodium hydroxide (NaOH) is prohibited. Use of corrosive products and/or instruments including abrasive sponges and metal brushes should be avoided. Cleaning must be performed by personnel trained in the general procedures involving contaminant removal. Automated washer/disinfectant systems are not recommended as the sole cleaning method for surgical instruments. An automated system may be used in addition to the following manual cleaning procedure.

1. Thoroughly clean all instruments and trays prior to use and as soon as possible after use (within a maximum of 2 hours post-operation) with intensive rinsing under cool tap water (<40°C) to remove gross soil. Do not allow blood and debris to dry on the instruments. If cleaning must be delayed, place instruments in a covered container with appropriate detergent (Enzol® Enzymatic Detergent or equivalent) to delay drying.
2. No instruments within this system require disassembly as part of the cleaning process.
3. The following table describes the required steps for thoroughly cleaning the system instruments.



Xpress™ Minimally Invasive Pedicle Screw System

Table 1: Cleaning Parameters

Step	Agent	Time (mm:ss)
	Instructions	
1. Initial Clean	Enzol Enzymatic Detergent Solution (or equivalent)	3:00
	<p>Add one (1) ounce (30 mL) of Enzol to one (1) gallon (3.8 L) of tap water. Soak instruments immediately after use and flush detergent through all channels until evidence of organic material is removed. Soak for three (3) minutes. While submerged, use a soft bristle brush that is intended for use with surgical instruments, to gently remove visible debris. Pay close attention to threads, crevices, lumens and hard to reach areas. If organic material is dried-on, extend soak time and use two (2) ounces (60 mL) of Enzol per one (1) gallon (3.8 L) of warm tap water.</p> <p>A cannula cleaner stylet is provided to be pushed thru cannulated instruments to clear out contaminants.</p>	
2. Rinse	Deionized water	1:00
	Thoroughly rinse each instrument with deionized water including all channels to remove detergent for one (1) minute.	
3. Ultrasonic Clean	Enzol Enzymatic Detergent Solution (or equivalent)	10:00
	Prepare a fresh solution by adding one (1) ounce (30 mL) of Enzol and one (1) gallon (3.8 L) of warm tap water to a sonication unit (Branson Bransonic® Ultrasonic Cleaner or equivalent). Fully immerse the instruments in the solution and sonicate for ten (10) minutes.	
4. Ultrasonic Rinse	Deionized water	1:00
	Thoroughly rinse each instrument with deionized water including all holes and cannulations to remove detergent for one (1) minute.	
5. Inspection	Unaided eye	
	Inspect each instrument for evidence of organic material. Repeat the ultrasonic clean and rinse steps if needed.	



Xpress™ Minimally Invasive Pedicle Screw System

4. If instruments are not dry, use filtered pressurized air or lint-free wipes to dry.

INSPECTION

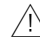
1. Carefully inspect each instrument and tray to ensure all visible blood and soil has been removed. Repeat cleaning if necessary.
2. Inspect instruments and instrument cases for damage. Check action of moving parts to ensure proper operation.
3. If damage or wear is noted that may compromise the proper function of the instrument or instrument case, do not use and contact customer service or your Xtant Medical representative for a replacement.
4. If corrosion is noted, do not use and contact customer service or your Xtant Medical representative for a replacement.

Manufacturer:



Xtant Medical

664 Cruiser Lane
Belgrade, MT 59714 USA
Phone: (888) 886-9354
Fax: (406) 388-3380

 **CAUTION:** Federal Law (USA) restricts these devices to use by or on the order of a physician.

REF

Part Number
X073-2003-EN Rev G
08/2020